

DICOM Conformance Statement

Web DICOM Viewer

Zargantana / Web DICOM Viewer


September 25, 2026

1. Introduction

This document is the DICOM Conformance Statement for the **Web DICOM Viewer**. This software component is designed for displaying DICOM images directly in a web browser. It supports uploading image files through a browser interface, visualizing them frame by frame with interactive navigation and rendering options.

2. Implementation Model

2.1 Application Data Flow

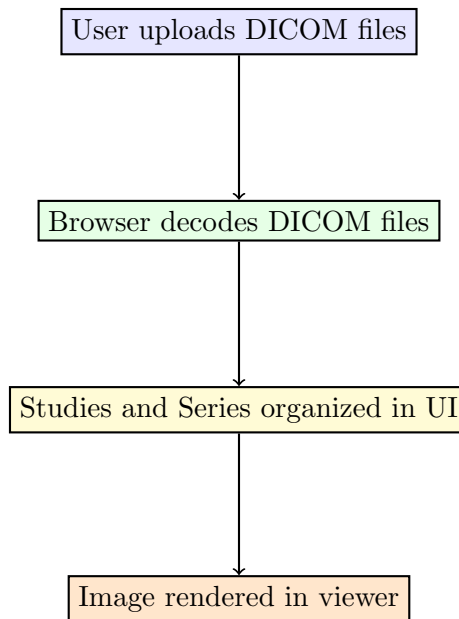

2.2 Functional Definition

The Web DICOM Viewer is a **VIEWER** application that does not establish network communication via DICOM protocol (no SCU/SCP roles). It operates solely on files provided through user upload.

3. AE Specifications

3.1 Web DICOM Viewer AE

- **SOP Classes Supported:**

- All storage SOP classes containing image data (e.g., CT, MR, XA, US, NM, CR, OT).
- **Supported Transfer Syntaxes** (decoding only; all verified with a synthetic regression suite):
 - Implicit VR Little Endian (1.2.840.10008.1.2)
 - Explicit VR Little Endian (1.2.840.10008.1.2.1)
 - Encapsulated Uncompressed Explicit VR Little Endian (1.2.840.10008.1.2.1.98)
 - Deflated Explicit VR Little Endian (1.2.840.10008.1.2.1.99)
 - Explicit VR Big Endian, retired (1.2.840.10008.1.2.2)
 - JPEG Baseline, Process 1 (1.2.840.10008.1.2.4.50)
 - JPEG Extended, Process 2 & 4, 8 and 12 bit (1.2.840.10008.1.2.4.51)
 - JPEG Extended, Process 3 & 5, arithmetic coding, 8 bit, retired (1.2.840.10008.1.2.4.52)
 - JPEG Spectral Selection, Non-Hierarchical, Process 6 & 8 and 7 & 9, 8 bit, retired (1.2.840.10008.1.2.4.53, .54)
 - JPEG Full Progression, Non-Hierarchical, Process 10 & 12 and 11 & 13, 8 bit, retired (1.2.840.10008.1.2.4.55, .56)
 - JPEG Lossless, Non-Hierarchical, Process 14, any predictor, 2–16 bit (1.2.840.10008.1.2.4.57)
 - JPEG Lossless, Non-Hierarchical, First-Order Prediction, Process 14 Selection Value 1 (1.2.840.10008.1.2.4.70)
 - JPEG-LS Lossless and Near-Lossless (1.2.840.10008.1.2.4.80, .81)
 - JPEG 2000 Image Compression, Lossless Only and lossy (1.2.840.10008.1.2.4.90, .91)
 - JPEG 2000 Part 2 Multi-component, Lossless Only and lossy (1.2.840.10008.1.2.4.92, .93), except array-based multi-component transformations (see Limitations)
 - High-Throughput JPEG 2000, Lossless Only, with RPCL Options (Lossless Only) and lossy (1.2.840.10008.1.2.4.201, .202, .203)
 - RLE Lossless (1.2.840.10008.1.2.5)
- **Supported Photometric Interpretations:** MONOCHROME1, MONOCHROME2, PALETTE COLOR (including segmented palette color lookup tables), RGB, YBR_FULL, YBR_FULL_422, YBR_PARTIAL_422, YBR_PARTIAL_420, YBR_ICT and YBR_RCT (JPEG 2000), and the retired HSV, ARGB and CMYK. Planar Configuration 0 and 1 are supported.
- **Image Display Capabilities:**
 - Grayscale pipeline: Rescale Slope/Intercept (linear Modality LUT), VOI (Window Center/Width with the LINEAR, LINEAR_EXACT and SIGMOID functions, or a VOI LUT), and presentation (MONOCHROME1 and Presentation LUT Shape INVERSE are inverted).
 - When the file defines several windows, the user can select any of them; the first one is used by default. Without a window or VOI LUT, images of up to 8 bits are displayed over their full range and deeper images use the minimum and maximum of the frame.
 - Pixel Padding Value / Range Limit: padding is shown in black and excluded from the automatic window.

- Supplemental Palette Color Lookup Tables (MONOCHROME2 with Pixel Presentation COLOR or MIXED).
- Multi-frame images are displayed frame by frame and played as a cine loop using Frame Time.

3.2 Supported Modalities

Modality Code	Description
CR	Computed Radiography
US	Ultrasound
XA	X-Ray Angiography
NM	Nuclear Medicine (e.g., gammagraphy)
OT	Other (Generic / Non-standard)
DX	Digital Radiography
CT	Computed Tomography (if uploaded)
MR	Magnetic Resonance (if uploaded)

4. Networking

Note: This application does not support DICOM network communication (no DICOM Upper Layer Protocol, DICOMweb, or DIMSE services).

5. Browser Interface and User Interaction

- Users upload DICOM files via an HTML file input element.
- The system displays a list of studies and associated series.
- Users can select a study and navigate through its images.
- Two display modes, chosen per series:
 - CT, MR and PT series with more than 10 images, and series of other modalities with more than 25 images: stack navigation with the mouse wheel.
 - Otherwise (always for US, CR, DX and XA): an image list with selection and, for multi-frame images, cine loop playback.

6. Limitations

- Not supported: JPEG 2000 Part 2 array-based multi-component transformations, hierarchical JPEG (processes 16–29), 12-bit JPEG with arithmetic coding or progression, JPEG XL, MPEG-2, MPEG-4 AVC/H.264 and HEVC/H.265 video, JPIP referenced pixel data, Float and Double Float Pixel Data, and ICC profile color management.
- JPEG Lossless with arithmetic coding (Process 15, 1.2.840.10008.1.2.4.58) is decoded on a best-effort basis for 8-bit images and has not been verified.
- Non-linear Modality LUT Sequences and per-frame functional group windows of enhanced multi-frame images are not applied.
- No support for DICOM C-FIND, C-MOVE, or other DICOM services.
- No integration with PACS or DICOMweb servers.

- Only file-based workflows supported.

7. Security

The viewer operates entirely on the client side, without uploading data to a remote server. Security and confidentiality of patient data depend on user handling and browser environment.

8. IHE Integration Statement

This viewer does not currently support IHE Integration Profiles such as RAD-68 or XDS-I.b.

9. Support and Contact Information

Web DICOM Viewer

Website: <https://visordicom.es>

Email: zargantana.reader@gmail.com